



PCT/GB 2004 / 0 0 2 7 6 4



INVESTOR IN PEOPLE

The Patent Office

Concept House

Cardiff Road

Newport

South Wales

NP10 8QQ

REC'D 21 JUL 2004

WIPO

PCT

## PRIORITY DOCUMENT


SUBMITTED OR TRANSMITTED IN  
COMPLIANCE WITH RULE 17.1(a) OR (b)

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

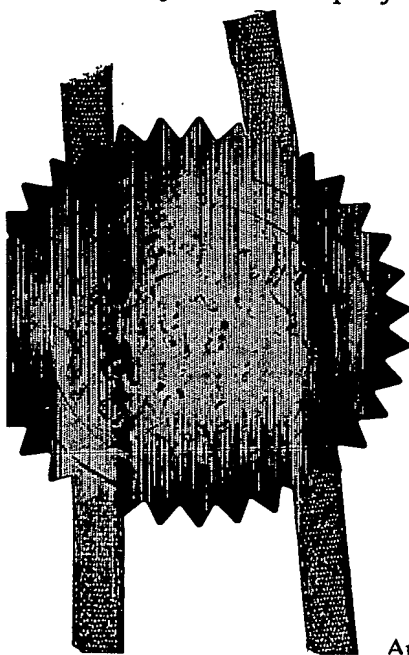

Signed

*Andrew Geney*

Dated 14 July 2004

The  
Patent  
Office

1/77

The Patent Office  
Cardiff Road  
Newport  
South Wales  
NP10 8QQ

**Request for grant of a patent**

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)



1. Your reference 3.81539/001

2. Patent application number  
(The Patent Office will fill in this part) 0328397.5

- 8 DEC 2003

3. Full name, address and postcode of the  
or of each applicant (underline all surnames)  
Dynal Biotech ASA  
Postboks 114 Smestad  
N-0309 Oslo  
Norway

Patents ADP number (if you know it)

816 16 55001

If the applicant is a corporate body, give  
country/state of its incorporation

Norway

4. Title of the invention Particles

5. Name of your agent (if you have one) Frank B. Dehn & Co.

"Address for service" in the United Kingdom  
to which all correspondence should be sent  
(including the postcode)

179 Queen Victoria Street  
London  
EC4V 4EL

Patents ADP number (if you know it)

166001 ✓

6. Priority: Complete this section if you are declaring priority from one or more earlier patent applications, filed in the last 12 months	Country	Priority application number (if you know it)	Date of filing (day / month / year)

7. Divisionals, etc: Complete this section only if this application is a divisional application or resulted from an entitlement dispute (see note f)	Number of earlier UK application	Date of filing (day / month / year)

8. Is a Patents Form 7/77 (Statement of  
inventorship and of right to grant of a patent)  
required in support of this request?

Answer YES if:

yes

- a) any applicant named in part 3 is not an inventor, or
- b) there is an inventor who is not named as an  
applicant, or
- c) any named applicant is a corporate body.

Otherwise answer NO (See note d)

# Patents Form 1/77

9. Accompanying documents: A patent application must include a description of the invention. Not counting duplicates, please enter the number of pages of each item accompanying this form:

Continuation sheets of this form

Description 16

Claim(s) 2

Abstract -

Drawing(s) -

10. If you are also filing any of the following, state how many against each item.

Priority documents -

Translations of priority documents -

Statement of inventorship and right to grant of a patent (Patents Form 7/77) -

Request for a preliminary examination and search (Patents Form 9/77) -

Request for a substantive examination (Patents Form 10/77) -

Any other documents (please specify) -

11.

I/We request the grant of a patent on the basis of this application.

Signature

*Neil Campbell*

Date 8 December 2003

12. Name, daytime telephone number and e-mail address, if any, of person to contact in the United Kingdom

Neil Campbell  
020 7206 0600

## Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

## Notes

- If you need help to fill in this form or you have any questions, please contact the Patent Office on 0645 500505.
- Write your answers in capital letters using black ink or you may type them.
- If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s) of the form. Any continuation sheet should be attached to this form.
- If you have answered 'Yes', Patents Form 7/77 will need to be filed.
- Once you have filled in the form you must remember to sign and date it.
- Part 7 should only be completed when a divisional application is being made under section 15(4), or when an application is being made under section 8(3), 12(6) or 37(4) following an entitlement dispute. By completing part 7, you are requesting that this application takes the same filing date as an earlier UK application. If you want the new application to have the same priority date(s) as the earlier UK application, you should also complete part 6 with the priority details.

## Improvements in Magnetic Polymer Particles

This invention relates to magnetic polymer particles carrying a chelating matrix loaded with a metal as well as to a process for the preparation of magnetic polymer particles carrying said chelating matrix. In particular, the invention relates to magnetic polymer particles carrying a carboxymethylated asparate (Cm-Asp) chelating group and to the coupling of the Cm-Asp group with the magnetic polymer particle.

Magnetic polymer particles are of general utility in various medical and biochemical fields, for example as transport vehicles for the delivery of pharmaceutical products, for diagnostic purposes, for separation and for synthetic purposes. Such particles rely upon their magnetic properties in order to perform these functions: in diagnostic assay applications, for example, application of a magnetic field to a sample containing an analyte bound to magnetic polymer particles allows the isolation of the analyte without the use of centrifugation or filtration; and in therapeutic applications, for example, application of a magnetic field to the patient may serve to target drug-carrying magnetic polymer particles to a desired body site.

By magnetic is meant herein that the polymer particles contain superparamagnetic crystals. Thus the magnetic polymer particles are magnetically displaceable but are not permanently magnetizable. Many processes for preparing magnetic polymer particles are known, a large number of which involve preparing maghemite- or magnetite-containing polymer particles from pre-formed magnetic iron oxides, e.g. magnetite. Some of processes involved are described in US-A-4,654,267 (Ugelstad) the contents of which are incorporated herein by reference.

The use of immobilised metal ion affinity chromatography (IMAC) has been known for many years. The

IMAC purification process is based upon the employment of a chelating matrix loaded with transition metal ions such as  $\text{Cu}^{2+}$  or  $\text{Ni}^{2+}$  which is capable of binding electron donating groups present on the surface of proteins, in particular the imidazole side chain of histidine. The electron donating group is believed to coordinate to vacant coordination sites around the metal ion. The interaction between the metal ion and the electron donating groups present on the protein surfaces can be altered by, for example, varying pH and hence purification can be achieved via reversible metal complex/protein interaction. Most commonly, if a protein is bound to a solid phase via the interaction between the metal ion and the imidazolyl side chain of histidine, the protein can be removed by addition of imidazole itself, i.e. by competitive elution.

Several different chelating ligands have been employed in IMAC to purify proteins. Nitrilo triacetate (NTA) (a tetradentate ligand) and the pentadentate ligand tris(carboxymethyl)ethylenediamine are examples of such ligands but these suffer from various disadvantages such as unspecific protein interaction, metal leakage etc.

US 6242581 proposes a solution to the metal leakage problem by the use of a carboxymethylated aspartate (Cm-Asp) group in IMAC where the bound transition metal ion has octahedral geometry. The ligand is said to be ideal for isolating histidine tagged recombinant proteins. Other advantages of Cm-Asp are discussed in US 5962641, e.g. resistance to reducing agents.

In these Patents the Cm-Asp ligand is bound to an agarose solid phase which is preferably cross-linked although other polymer matrices such as polystyrene, nylon and SEPHAROSE are suggested. Whilst these matrices may be magnetic the magnetic particles do not remain in suspension and the solid phases are therefore of limited use in assays.

It has now been surprisingly found that the Cm-Asp


chelating ligand can be coupled to a magnetic polymer particle giving rise to a moiety that possesses not only the ability to bind histidine-tags in recombinant proteins but also magnetism thereby allowing the skilled biochemist more flexibility in his assaying procedures. The inventors have also devised ways to couple the Cm-Asp ligand to the magnetic polymer particles in high yield thereby producing an excellent IMAC agent.

Viewed from a first aspect, therefore, the present invention provides a magnetic polymer particle bound to a carboxymethylated aspartate chelating ligand.

Viewed from another aspect the invention provides a magnetic polymer particle bound to a carboxymethylated aspartate ligand chelating a metal ion.

Viewed from another aspect the invention relates to a process for the preparation of a magnetic polymer particle as hereinbefore defined comprising reacting a magnetic polymer particle with a Cm-Asp chelating ligand.

The Cm-Asp ligand bound to the magnetic polymer particle (MPP) is depicted below both in its uncoordinated state and coordinated to a metal ion (the wavy line representing a bond or a linker between the Cm-Asp and particle):

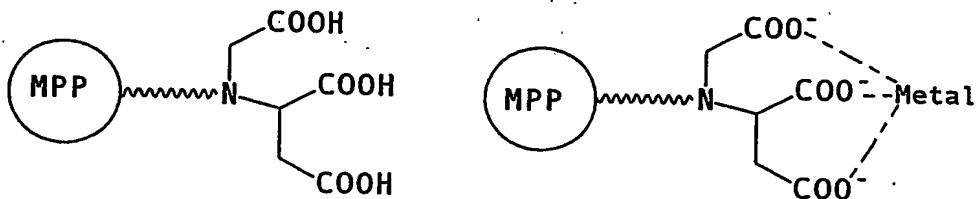

25

The magnetic polymer particles used in the process of the invention may be any magnetic polymer particle e.g. as described in US-A-4,654,267. The particles are preferably porous to allow the presence of the superparamagnetic crystals in the pores thereof. The surface of the magnetic

30

particles is normally functionalised to allow coupling of the Cm-Asp ligand to the polymer particle, e.g. it may be functionalised to carry any known surface structure such as carboxyl groups, tosyl groups, amino groups, epoxy groups, maleamido groups, thiol groups etc. Most preferably the surface is amine functionalized before Cm-Asp coupling. Alternatively, an amine functionalised surface can itself be further functionalised to carry other functional groups, e.g. COOH groups.

- 10       The polymer particle is preferably made from combinations of vinylic polymers (e.g. styrene), acrylates and/or methacrylates. The polymeric material may optionally be crosslinked, for example by incorporation of cross-linking agents, for example as comonomers, e.g. 15 divinylbenzene (DVB) or ethyleneglycol dimethacrylate. Appropriate quantities of the cross-linking agents (e.g. comonomers) required will be well known to the skilled man. Preferably the polymer is a cross-linked styrenic polymer (e.g. a styrene-divinylbenzene polymer, surface 20 functionalized by the use of a nitro-group containing comonomer, e.g. nitro-styrene, and subsequent reduction) or a cross-linked (meth)acrylic polymer surface functionalized by the use of an epoxy-group containing comonomer (e.g. glycidylmethacrylate) and subsequent amination (e.g. by 25 reaction with ethylene diamine).

- The superparamagnetic crystals in the polymer particles used in the process of the invention may be of any material capable of being deposited in superparamagnetic crystalline form in the porous polymer 30 particles. Magnetic iron oxides, e.g. magnetite or maghemite are preferred; however the crystals may be of mixed metal oxides or other magnetic material if desired. The total quantity of crystalline magnetic material present is generally more than 1%, preferably more than 3%, 35 desirably more than or equal to 5% (by weight, e.g. up to 40% wt. The percentage is calculated on a Fe (or
-

equivalent metal in the case of magnetic materials other than iron oxides) weight basis based upon the overall dry weight of the coated particles.

Polymer particles according to the various aspects of the present invention will generally have sizes (i.e. diameters) that are generally in the micrometer range, e.g. 0.3 to 100  $\mu\text{m}$ , especially 0.5 to 50  $\mu\text{m}$ , more especially 0.8 to 5  $\mu\text{m}$ , e.g. 0.8 to 1.2  $\mu\text{m}$ .

Typically the porous particles used will have a surface area of at least 15  $\text{m}^2/\text{g}$  (measured by the BET nitrogen absorption method), and more preferably at least 30  $\text{m}^2/\text{g}$ , e.g. up to 700  $\text{m}^2/\text{g}$ , when corrected to a mean particle diameter of 2.7  $\mu\text{m}$  (i.e. multiply surface area by 2.7/MD, where MD is the mean diameter in micrometers). Similarly scaled, the particle pore volume is preferably at least 0.1 mL/g.

Typically, the polymer particles are spherical and substantially monodisperse before they are coated and especially preferably remain spherical and substantially monodisperse once they have been coated.

By substantially monodisperse it is meant that for a plurality of particles (e.g. at least 100, more preferably at least 1000) the particles have a coefficient of variation (CV) of less than 20%, for example less than 15%, preferably less than 12%, more preferably less than 11%, still more preferably less than 10% and most preferably no more than about 8%, e.g. 2 to 5%. CV is determined in percentage as

$$\text{CV} = \frac{100 \times \text{standard deviation}}{\text{mean}}$$

where mean is the mean particle diameter and standard deviation is the standard deviation in particle size. CV is preferably calculated on the main mode, i.e. by fitting a monomodal distribution curve to the detected particle size distribution. Thus some particles below or above mode size



may be discounted in the calculation which may for example be based on about 90% of total particle number (of detectable particles that is). Such a determination of CV is performable on a Coulter LS 130 particle size analyzer.

5        Functionalisation of the polymeric material may take place after polymerisation by, for example, nitration and subsequent reduction of the thus-formed nitro groups to pendant amine groups; or direct amination, for example by treatment with amino ethanol. As further alternatives,  
10       polymeric particles prepared by the well-known Ugelstad two-step swelling process and the improvements thereto disclosed in WO 00/61647 (Dyno) may be used. Porous polymer particles produced according to the processes described in this publication may have magnetic particles  
15       deposited in their pores by standard techniques.

As a further possibility, porous polymer particles may be prepared from nitro styrene and DVB, and magnetic material introduced as taught in US-A-4,654,267.


20       The superparamagnetic polymer beads sold by Dynal Biotech ASA under the trade names Dynabeads, especially Dynabeads MyOne are especially preferred. Dynabeads are particularly advantageous since they remain in suspension and do not exhibit magnetic particle sedimentation often associated with other magnetic beads. Dynabeads also show  
25       excellent magnetic mobility compared to other magnetic particles in which high levels of iron are present. Dynabeads exhibit beneficial kinetics allowing shorter reaction times and higher throughputs. Their unspecified binding is lower than other magnetic beads and their proper  
30       use results in a concentration of the desired material taking place resulting in easier and more efficient washing procedures. Finally Dynabeads, e.g. MyOne beads are easy to automate and are monodisperse.

35       The Cm-Asp ligand is bound to the magnetic polymer particle. By bound is meant that the ligand is covalently linked to the polymer particle, optionally using a linking

---

group as discussed in detail below. The Cm-Asp ligand can be bound to the magnetic polymer particle by various procedures although it is preferred if there are at least three linking atoms between the polymer particle surface and the nitrogen atom of the Cm-Asp, e.g. the styrene surface and the nitrogen atom of the Cm-Asp ligand. Preferably there are at least 6 atoms separating the Cm-Asp ligand from the magnetic polymer particle surface, more preferably there are between 6 and 20 atoms separating the Cm-Asp ligand from the magnetic polymer particle surface.

Hence, in a preferred embodiment the invention provides a composition of formula (I)

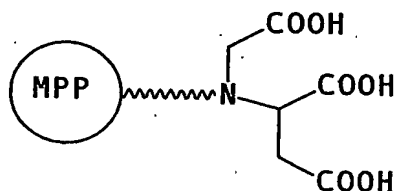

15

(where MPP is a magnetic polymer particle and the wavy line represents a linking group comprising at least three atoms) or an analogue thereof in which a metal ion is chelated.

In US 6242581 aspartic acid is coupled to the solid phase prior to carboxymethylation to form the Cm-Asp ligand however it has not been possible to use this technique to provide a Cm-Asp group on a magnetic polymer particle. Rather, the inventors have devised alternative syntheses in which the Cm-Asp ligand is fully formed prior to coupling to the magnetic polymer particle.

In this regard, it has been found that when there are fewer than 3 atoms between the polymer surface and Cm-Asp ligand then coupling yields are low. In contrast to an agarose support carrying Cm-Asp (as describe in US-A-5962641), it is necessary in the present invention to ensure that coupling yields between the magnetic polymer particle and Cm-Asp are relatively high. The surface area

of an agarose support is considerably greater than that of a polymer particle and hence the binding of Cm-Asp to the support does not need to be achieved in high yield for a useful IMAC chelating agent to result. In the present case, yields need to be much higher to ensure that enough polymer particles carry the Cm-Asp ligand and hence to ensure that IMAC can be successfully carried out.

It is preferred if the at least 3 atom linker comprises an amino group (-NH-). Magnetic polymer beads are often made from styrene polymers which are nitrated to form NO<sub>2</sub> groups on the surface. After reduction of these groups, e.g. using ammonia, amino groups are formed and these form the most common link from the polymer particle surface.

The next portion of the linker preferably represents the residue of an electrophile; i.e. the group which remains after reaction of the electrophile with a nucleophile. Hence, the linker may comprise an oxo group (C=O, the residue of an ester/carboxyl group), a -CH(OH)CH<sub>2</sub>- group (the residue of an epoxide), -CH<sub>2</sub>- (where the electrophile is, for example a CH<sub>2</sub>Hal). The linker may also incorporate a number of atoms linking the actual electrophile to the -NH- group, e.g. an alkylene chain or ether chain, e.g. as in -CH<sub>2</sub>CH<sub>2</sub>CH<sub>2</sub>-, or -CH<sub>2</sub>CH<sub>2</sub>CH<sub>2</sub>-O-.

A final portion of the linker represents the residue of a nucleophile from the Cm-Asp, i.e. the residue which results after reaction of this nucleophile with the electrophile. As discussed in more detail below this may be a aminoalkylene or aminoether/polyether, thiol or hydroxyl residue.

Hence the wavy line in formula (I) can represent -NH-L<sub>1</sub>-Er-Nr-L<sub>2</sub>- wherein L<sub>1</sub> represents a 1 to 10 atom linker to the electrophile residue (Er), and L<sub>2</sub> represents a 1 to 10 atom linker to the nucleophile residue (Nr).

It is of course within the scope of the invention for the magnetic polymer particle to carry a nucleophile with

the Cm-Asp being functionalised to carry an electrophilic group.

In a preferred embodiment the polymer particle should be functionalised to carry a coating which can react with the Cm-Asp ligand to couple the magnetic particle to the Cm-Asp.

In an especially preferred embodiment, a particle coating is provided which carries a carbon-carbon double bond. This can be achieved by, for example, reaction of the particle with an allyl or vinyl compound, e.g. butenoic acid. Hydroxy functionalised particle surfaces can be reacted with allyl bromide to form double bonds on the particle surface. Also, carboxy functionalised particle surfaces can be reacted with allyamines to provide double bonds on the particle surface. The Cm-Asp may then be coupled directly to the double bond using appropriate chemistry or more preferably, the double bond may then be reduced e.g. in the presence of aqueous halide to provide a halide electrophile which can be reacted with the Cm-Asp ligand to ensure successful coupling.


Another preferred preparation process involves functionalising the surface of the magnetic polymer particle to carry carboxyl groups. The carboxylic acid groups can be activated by reaction with N-hydroxysuccinimide esters and reacted with a Cm-Asp ligand as discussed above.

The Cm-Asp ligand can coordinate any metal ion. By metal ion is meant any ion of a metal from groups 1 to 13 of the periodic table, a lanthanide or actinide ion or an ion of Si, Ge, Sn, Pb, As, Sb, Bi, Te, Po or At. The metal ion should preferably be a transition metal ion and should preferably be in the 2+ oxidation state. Preferred metals are Ni, Fe, Ga, Mn, Co, Cu and Zn of which Co, especially  $\text{Co}^{2+}$  is more preferred. Coordination can be easily effected by exposing the Cm-Asp to, for example, the metal (II) chloride.

The Cm-Asp ligand may too be functionalised prior to coupling with the magnetic polymer particle. For example, it has proved advantageous to provide the Cm-Asp ligand with a linking group carrying a primary nucleophile to aid reaction with electrophilic groups on the particle surface. The nitrogen atom of the Cm-Asp ligand is secondary and it has been found that this atom is too unreactive, perhaps due to steric hindrance, to react in high yield with electrophilic groups, e.g. halides, on the particle surface.

It is preferred therefore to couple the Cm-Asp to a linker group having at least two atoms and comprising a nucleophile such as an amine, hydroxyl or thiol group. Preferably the linker is an alkylamine, e.g. C<sub>5</sub>/6-alkylamine linker or an ether/polyether linkage e.g. comprising one or two oxygen atoms and 3 to 6 carbon atoms. Coupling of the linker to the Cm-Asp (via the nitrogen atom thereof) is achieved using known chemistry as described in the Examples.

Thus viewed from another aspect the invention provides a process for the preparation of a magnetic polymer particle bound to a Cm-Asp ligand comprising reacting a Cm-Asp ligand of formula (II)

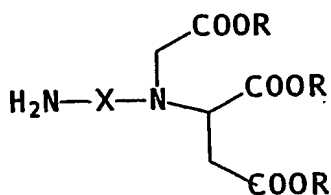


25

(wherein each R independently represents hydrogen or a protecting group and X represents a 2 to 10 atom linker, e.g. an C<sub>2-10</sub> alkylene linker, especially a C<sub>5/6</sub>-alkylene linker with a magnetic polymer particle functionalised to carry an electrophilic coating, e.g. an ester, epoxide, allyl, alkyl halide etc coating. Compounds of formula (II)

30

are themselves new and form a further aspect of the invention along with the Cm-Asp ligand itself, i.e. a compound of formula (III)

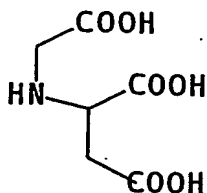

5

In some embodiments of the invention it may be necessary to protect the carboxyl groups of the Cm-Asp ligand during syntheses. This can be easily effected using known protection strategies, e.g. using an ester protecting group which can be hydrolysed in acid or base as is known in the art.

The magnetic polymer particles carrying the Cm-Asp ligand with associated metal ion can in general be used for attaching to and combining with peptides, proteins or other polymers (e.g. antibodies) and are hence of use in a wide variety of assays. They are of particular use, however, in the isolation of His-tags in recombinant proteins. Hence viewed from another aspect the invention provides the use of a magnetic polymer particle bound to a Cm-Asp ligand, said ligand coordinating a transition metal ion, in an assay. Suitable assays and ways to carry these out are known by the skilled biochemist.

For example, the capture of histidine-tagged proteins on the Cm-Asp functionalised particles of the invention has various applications. The rapid reaction kinetics and gentle handling of isolated proteins make this technology well suited for the "pull down" of large protein complexes. Thus Cm-Asp functionalised beads may be used in sample preparation for mass spectrometry analysis. It is believed that complexes isolated with the Cm-Asp functionalised beads may be more intact than complexes isolated with

30

columns or other solid supports including other magnetic particles with uneven surfaces and are therefore ideal for use in mass spectrometry sample isolation.

5 The Cm-Asp technology may also act as a solid phase for use in assay procedures. The Cm-Asp beads are not prone to aggregation and are highly dispersed in solution and show a low degree of non-specific binding. These properties allow for high quality screening results and protocols that are easily automated on a wide range of  
10 automation platforms. The Cm-Asp beads may also be used in phage display perhaps as a solid phase or to purify expressed phage display selected proteins from a library.


In general therefore the capture of histidine tagged proteins may allow microscale protein purification, clean  
15 up of mutated protein libraries, denaturing elution of protein/peptide, mild elution of proteins/peptide, protein-protein interaction studies and screening technologies, e.g. for drug discovery, molecular display, aptamer screening, phage display, engineered enzyme screening and  
20 diagnostics.

The invention will now be described further by reference to the following non-limiting examples.

---

# Reactant Preparation

The Cm-Asp triester below is prepared as follows:

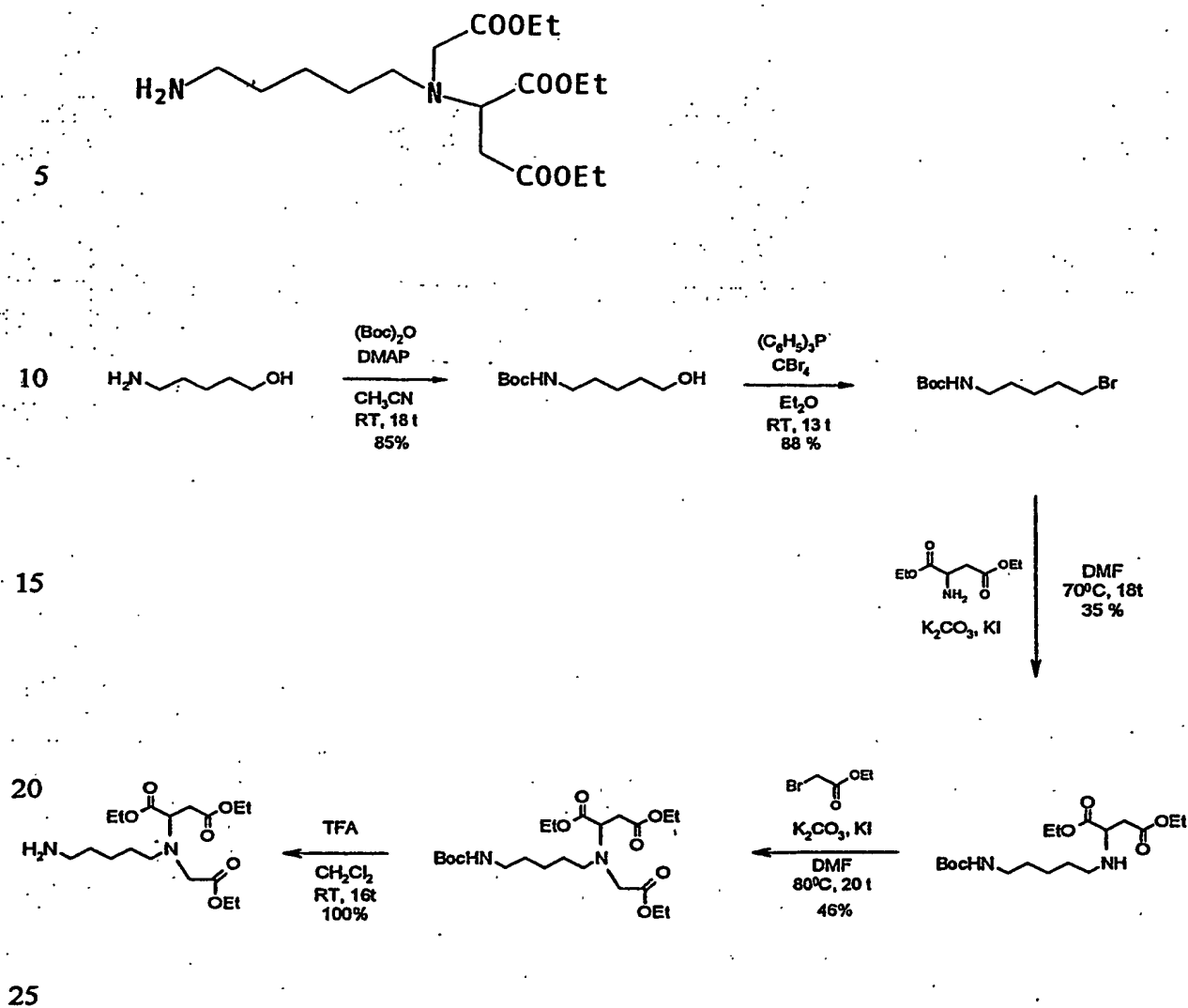

## Example 1:

### Bromination

17.3 g of a methanol suspension of the magnetic styrene particles having 0.5 mmol/g allyl groups are washed four times with 45 mL sodium acetate buffer (pH = 5.9). After



adjusting the particle content to 9 wt%, 0.96 g of pyridinium tribromide dissolved in 10 mL DMF is added while stirring at 350 rpm. After five minutes at room temperature the particles are washed five times with 45 mL deionised water.

Example 2:

Functionalization with Cm-Asp chelator

18.0 g of a suspension of the particles prepared as in Example 1 are washed three times with 20 mL of 50mM sodium bicarbonate. The particle content is adjusted to 12 wt%. To the suspension 0.17 g of the Cm-Asp triester (prepared as described above) is added. 50mM sodium bicarbonate is added until a particle content of 10 wt% was achieved. The reaction mixture is shaken at 600 rpm at 40°C for 15 hours. The particles are then washed four times with 20 mL deionised water.

Example 3:

Hydrolysis

20.0 g of a suspension of particles prepared as in Example 2 are washed twice with 20 mL of 1M lithium hydroxide. After adjusting the particle content to 10 wt% the mixture is shaken at 250 rpm for four hours at room temperature. The particles were then washed with deionised water until pH 6-7.

Example 4:

Cobalt-loading

250 mg of particles prepared as in Example 3 are washed twice with 5 ml reverse osmosis-water. 5 ml 2,5 mM  $\text{CoCl}_2$  are added to the particles and incubated for 5 h. The tube

---

is placed in a magnet, and the supernatant is removed. The particles are washed twice with 5 ml phosphate buffered saline (0,01% Tween 20, pH 7,4). The particles are then washed once in 20% ethanol. The particles are stored in 20% ethanol.

#### Example 5

##### Functionalisation of carboxylic acid groups to N-hydroxysuccinimide ester

10

50 g of a suspension of 5.0 g of the particles of MyOne Carboxylic acid beads are acidified by washing with 0.1 M acetic acid (3 x 50 mL). The acidified particles (which have a carboxylic acid content of 0.5 mmole/g DS) are then washed with acetone (4 x 50 mL) and concentrated on a magnet. Extra acetone is added until a total of 35.6 g suspension is achieved. N-hydroxysuccinimide (2.90 g, 25 mmole) and diisopropylcarbodiimide (3.16 g, 25 mmole) are then added. The reaction mixture is stirred at room temperature for 5 hours. The particles are then washed with acetone (5 x 50 mL).

15

20

#### Example 6:

##### Functionalization with Cm-Asp chelator

25

44 g of an acetone suspension of the beads of Example 5, are washed three times with 50 mL isopropanol. After adjusting the particle content to 12 wt%, 5.6 g of triethylamine is added. 0.10 g of the Cm-Asp triester (prepared as described above) dissolved in isopropanol, is then added. This results in a particle content of 10 wt%. The reaction mixture is then shaken at 250 rpm at room temperature for 20 hours. The particles are washed three times with 50 mL of isopropanol.

30

35

**Example 7:****Functionalization with Cm-Asp chelator and ethanolamine**


To 10 g of an isopropanol suspension of the particles  
5 prepared as in Example 6, 0.32 g of ethanolamine is added.  
The reaction mixture is then shaken at 250 rpm at room  
temperature for 18 hours. The particles are then washed  
three times with 10 mL of isopropanol.

**10 Example 8****Functionalization with Cm-Asp chelator**

1,2 gram of dry Dynabeads 270 Epoxy are mixed with 8,8 gram  
of 50 mM sodium bicarbonate. To the suspension 0,17 grams  
15 of the Cm-ASP triester (prepared as described above) are  
added, and the reaction mixture is shaken at 600 rpm at  
60°C for 16 hours. The particles are worked up by washing  
four times with 20 ml deionised water.


## Claims

1. A magnetic polymer particle bound to a carboxymethylated aspartate chelating ligand.
- 5 2. A magnetic polymer particle bound to a carboxymethylated aspartate ligand chelating a metal ion.
- 10 3. A process for the preparation of a magnetic polymer particle as hereinbefore defined comprising reacting a magnetic polymer particle with a Cm-Asp chelating ligand.
- 15 4. A process for the preparation of a magnetic polymer particle bound to a Cm-Asp ligand comprising reacting a Cm-Asp ligand of formula (II)

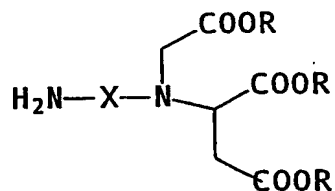

- 20 (wherein each R independently represents hydrogen or a protecting group and X represents a 2 to 10 atom linker) with a magnetic polymer particle functionalised to carry an electrophilic coating.


5. A compound of formula (II)

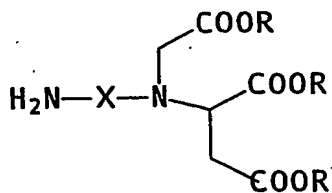

25

- (wherein each R independently represents hydrogen or a protecting group and X represents a 2 to 10 atom linker).

6. A compound of formula

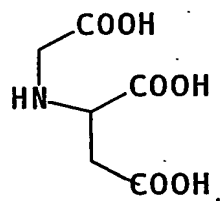

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☒ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☒ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**